


‘ Loughborough ‘ Glass Water Still

This traditional glass water still has a long history of reliability and robustness and has served countless laboratories worldwide for many years.

Incorporating a heavy duty chromium plated heating element with thermostatic protection the Still is supplied complete with bench or wall mountable stand. The robust borosilicate glass system features a boiling flask , combined condenser and vapour trap and water constant level device.

A boiler cleaning syphon tube enables the boiler contents to be drained to allow descaling so avoiding the need to dismantle the glassware.

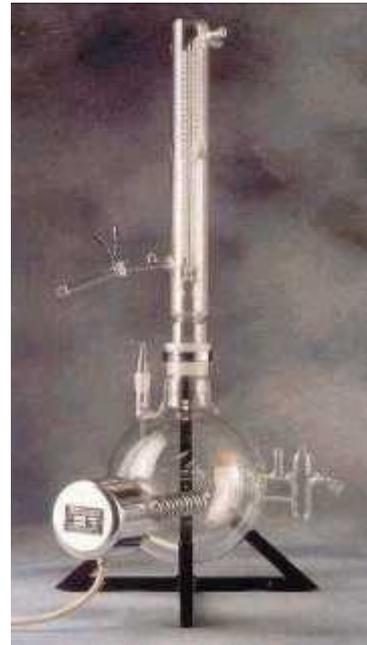

Output sizes of 4 and 8 lt/hr are offered. Supplied complete with all connection hoses and wall / bench mounting bracket.

Specifications & Ordering Information			
Catalogue No	Output	Power	Dimensions WxDxH (mm)
WGS-560-020M	4 L/hr Single	220V 3kW	450 x 390 x 820
WGS-560-040G	4 L/hr Single	240V 3kW	450 x 390 x 820
WGS-570-020F	8 L/hr Single	220V 6kW	580 x 430 x 1040
WGS-570-040W	8 L/hr Single	240V 6kW	580 x 430 x 1040

Monarch Way
Belton Park
Loughborough
LE11 5XG UK
Tel : +44 (0) 1509 224613 / 224615 Fax : +44 (0) 1509 260210
Email : sales@fistreem.co.uk
www.fistreeminternational .com



Fistreem International Ltd