
- High-temperature shaker—up to 85°C—available
- Assortment of platforms available for increased versatility
- Lifetime warranty on drive mechanism



The Lab-Line 3527 Series Bench Top Environmental Shakers offer the convenience of a bench top incubator combined with an orbital shaking mechanism. They are ideal for cell and bacterial culture studies, genetics research, and increasing solubility rates.

- ♦ Compact design conserves valuable bench top space
- ♦ Models 3527, 3527-1, 3528-5 and 3528-6 feature a temperature range from slightly above ambient to 65°C; Models 3527-5 and 3527-6 offer higher temperature range to 85°C.
- ♦ Chamber accommodates a variety of platforms, holding flasks up to 1<</p>
- ♦ Models 3528-5 and 3528-6 hold flasks up to 2<
- ↑ Triple eccentric drive powered by 1/15 H P motor rotates shaker platform through 3/4" (2 cm) diameter orbit
- ♦ Solid state controller varies speeds from 40 to 400 R P M
- ♦ Easy-to-read LED temperature display
- $\mbox{\Large Φ}$ Safety interlock halts shaker action when cover is opened



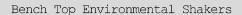

See page 104 for L-C Incubators

- ♦ Independent hydraulic thermostats maintain temperature with ±0.2°C control; provide automatic over-temperature protection
- ♦ Smooth counterbalanced drive is isolated from shaker chamber to prevent damage from spills or glassware breakage
- ♦ Drive rotates on nine permanently lubricated ball bearings

ACCESSORIES						
Model	Description					
CC*	Cooling Coil (factory installed)					
GM*	Gassing Manifold (factory installed)					
4628-40	Light Bank, 120V, 60 Hz					
4628-41	Light Bank, 240V, 50 Hz					
Note: If a 3527 with a cooling coil is required, order 3527CC.						

SPECIFICATIONS AND ORDERING INFORMATION															
Model	Volts*	Amps	Watts	Temp. Range	Temp. Control at 37°C (in flask)	Temp. Uniformity at 37°C (in flask)	Chamber Dimensions W" x D" x H" (cm)	Overall Dimensions** W" x D" x H" (cm)	Net Wt. Ibs. (kg)	Ship Wt. Ibs. (kg)					
3527	120	5.0	600	600	amb.+5°										
3527-1	240	2.5		to 65°C	±0.2°C	±0.5°C	18 x 18 x 10	201/8 x 27 x 191/4	110	150					
3527-5	120	8.3	1000	1000	1000	1000	1000	1000	amb.+5°	±0.2 C	10.5 C	$(46 \times 46 \times 25)$	(51 x 69 x 49)	(50)	(68)
3527-6	240	4.2		to 90°C											
3528-5	120	5.0	600	600	600	600	600	amb.+5°	.0.000	.0.0%	21 x 20 ³ / ₄ x 13 ³ / ₄	22 ³ / ₁₆ x 29 ³ / ₄ x 22 ⁵ / ₈	135	235	
3528-6	240	2.5		to 65°C	±0.2°C	±0.6°C	(53 x 53 x 35)	(56 x 76 x 56)	(61)	(107)					

*All models rated 50/60 Hz

^{**}Opened lid adds 16½" (42 cm) to height